



CoreLife

MEDICAL COMPRESSION STOCKINGS

Corelife fully meets the medical necessities for thigh high compression stockings. Besides being durability and robustness, the stockings are quite healthy, where such features are provided by state-of-the-art technology. Corelife ensures 100% that the applicable pressure is the right pressure, and the manufacturing is carried out in compliance with the standard Eko-Tex 100.

Entire thigh high stocking models feature silicon bands. The stockings are air permeable, antibacterial and hypoallergenic. The ones who suffer from diabetes or infection in feet, either for short or extended terms, may prefer open-toe models. Others may prefer closed toe models. Double weaving technology is employed for open-toe models. For open-toe models, an exclusive double layer weaving is utilized. Hence, neither incision nor unfavorable situation arises out on the feet or toes

EXPLANATION OF COMPRESSION CLASSES

Class 1 (CLL-1) – LOW PRESSURE COMPRESSION STOCKINGS – 18-21 mmHg

Class– 2 (CLL-2) – MEDIUM PRESSURE COMPRESSION STOCKINGS – 23-32 mmHg

Class –3 (CLL-3)– HIGH PRESSURE COMPRESSION STOKINGS - 34-46 mmHg

MATAREIL COMPOSITON

	Elastane (EL)	Polyamide (PO)
CLL-1	25 %	75 %
CLL-2	31 %	69 %
CLL-3	37 %	63 %

Under Knee Compression Stockings / Open Toe – Closed Toe



Above Knee Compression Stockings / Open Toe – Closed Toe

With silicone drop support



Pantyhose Compression Stockings / Open Toe – Closed Toe

With Belt



Medical Compression Stockings



Comfort For Legs





Ankle Measurement

1. Measure around your ankle (this measurement is called the circumference). Place the measuring tape at the narrowest part of the ankle. This must always be your first measurement.

Calf Measurement

2. Sit on a chair with your legs at a 90-degree angle. Find the widest part of your calf and measure the circumference.

Thigh Measurement

3. Find the widest part of your thigh - right under your buttocks. Measure the circumference of this part of your thigh.

Model	A-D Calf	A-G Thigh	A-T Panty Maternity Panty
Size	All Stockings Ankle Measurement (cm)	Calf Stockings Calf Measurement (cm)	Thigh Panty, Maternity Stockings Thigh Measurement (cm)
1	18 - 20	27 - 33	43 - 57
2	20 - 22	29 - 36	45 - 62
3	22 - 24	32 - 39	49 - 67
4	24 - 26	34 - 42	53 - 72
5	26 - 28	36 - 45	56 - 77
6	28 - 30	38 - 48	60 - 81
7	30 - 32	40 - 50	64 - 85

Proper compression treatment depends on proper measurement of ankle, calf and thigh.




Explanation Of Compression Classes	
CCL 1	Compression in mmHG 18.0 - 21.0
CCL 2	Compression in mmHG 23.0 - 32.0
CCL 3	Compression in mmHG 34.0 - 46.0

Material Composition		
	Elastane (EL)	Polyamide(PA)
CCL 1	25 %	75 %
CCL 2	31 %	69 %
CCL 3	37 %	63 %

Leg Circumference	SIZE (cm)						
	1	2	3	4	5	6	7
cG	43 - 48	45 - 52	49 - 56	53 - 60	56 - 64	60 - 68	64 - 72
cG	49 - 57	53 - 62	57 - 67	61 - 72	65 - 77	69 - 81	73 - 85
cG	43 - 57	45 - 62	49 - 67	53 - 72	56 - 77	60 - 81	64 - 85
cF	39 - 52	41 - 56	44 - 60	47 - 65	50 - 69	53 - 73	57 - 78
cE	30 - 37	33 - 40	35 - 43	37 - 45	39 - 48	41 - 51	42 - 53
cD	27 - 33	29 - 36	32 - 39	34 - 42	36 - 45	38 - 48	40 - 50
cC	28 - 34	30 - 37	33 - 40	35 - 43	37 - 46	39 - 49	41 - 51
cB1	23 - 27	24 - 29	26 - 32	29 - 35	31 - 37	33 - 39	35 - 41
cB	18 - 20	20 - 22	22 - 24	24 - 26	26 - 28	28 - 30	30 - 32

Length of Leg (cm)	
IG Normal	72 - 83
IG Short	62 - 71
ID Normal	39 - 44
ID Short	34 - 38

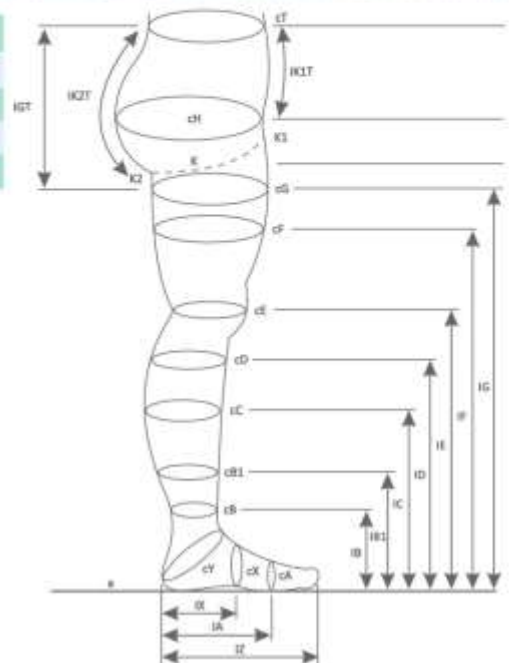

5747/11 Sok No 2 MTK Ssitesi
Bornova / Izmir /Turkey
+90 232 231 11 05
www.corelife.com.tr